



Roller and tilt mixer "Movil-Rod II"

FIXED SPEED 36 RPM.

FEATURES

Roll and tilt mixer with 5 rollers, double action, gently rolls the tubes while rocking, produces a sinusoidal movement, providing excellent sample homogenisation. Recommended for suspensions of biological and haematological fluids.

MODEL

Part No.	For tube Ø up to mm	Speed r.p.m.	Length (rollers) cm	Height / Width / Depth (exterior) cm	Power W	Weight Kg
7001726	35	36	23	14 26 15,5	6	3



NEW



Rocking Mixer "Vibromatic"

DIGITAL ELECTRONIC CONTROL OF SPEED AND TIME.



Optimum shaking frequency: 950 oscillations per minute

APPLICATIONS

Dissolution, extractions and homogenisation.

FEATURES

Electronic control from 100 to 950 oscillations per minute. (1900 vibrations) and timer from 1 to 999 minutes or continuous operation.

Oscillation amplitude: 8 mm.

Easily detachable arms equipped with 8 clamps that allow all types of tubes and flask with diameters of between 5 to 50 mm can be held in position.

Maximum load: 4 kg.


CONTROL PANEL

1. Digital display indicator of osc./min. and time.
2. Parameter display selector.
3. Parameter selector, speed or time.
4. Push button increase value.
5. Push button decrease value.
6. Push button start-stop.



MODEL

Part No.	Oscillations p.m.	Maximum capacity Kg	Height / Width / Depth (exterior) cm	Power W	Weight Kg
7000384	100 - 950	4	23 77 27	34	11

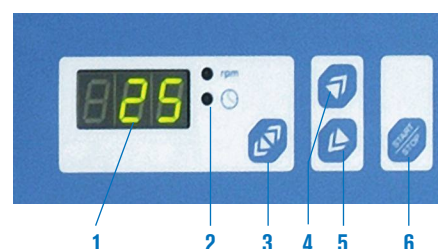

ACCESSORIES

2 in 1. Platform support for Erlenmeyer flasks.

Platform support for Erlenmeyer flasks, made from AISI 304 stainless steel with clip adapter position holes for 50 ml to 1000 ml flasks. Platform dimensions: 38 cm long x 19 cm wide.

Weight: 2 Kg.

Platform capacity:

- 18 Erlenmeyer flasks of 50 ml.
- 12 Erlenmeyer flasks of 100 ml.
- 8 Erlenmeyer flasks of 250 ml.
- 6 Erlenmeyer flasks of 500 ml.
- 3 Erlenmeyer flasks of 1000 ml.

Part No. **3001221**



Model Vibromatic with adapter platform.

Adapters for Erlenmeyer flasks.

Made from tempered stainless steel.

Part No.

- | | | |
|----------------|-------------|----------|
| 1001003 | Adapter for | 50 ml. |
| 1000978 | Adapter for | 100 ml. |
| 1000979 | Adapter for | 250 ml. |
| 1000980 | Adapter for | 500 ml. |
| 1000981 | Adapter for | 1000 ml. |

